



QUANTASEP® 1000LX

Dual Flow Range Scale-Up System

QuantaSep® 1000LX is a unique integrated system that enables low flow methods development and high flow scale up in one machine. With its powerful and easy to use Windows®-based software, it allows the user to automatically change buffers, collect fractions, and run gradients. The user can control each step in the purification based on UV, pH, conductivity, time, volume, air sensors or any external user-defined variable. All event, alarm, method and chromatography data is logged for analysis or archiving. When methods development is complete, simply switch to the higher flow tubing to scale-up and make clinical product.



Low Flow: up to 100 ml/min; High Flow: up to 1000 ml/min

8 Mission Critical Ways a QuantaSep® Can Help Get Your Department In Control


- ✓ Reduce Project Backlog
- ✓ Reduce Tedious Documentation
- ✓ Fast Development
- ✓ Easy Scale-Up 1000X!
- ✓ Achieve Clinical Deadlines
- ✓ Improve Process Economics by 50%
- ✓ Unattended Clinical Production 24/7 Now!
- ✓ Space Efficient

QuantaSep® 1000LX Highlights

- Automates buffer delivery, column switching, fraction collection, based on UV, pH, and conductivity.
- Small footprint can fit on bench-top, laminar flow hood or cold room.
- Dual dynamic flow range - 1-100ml/min & 10-1000 ml/min
- Graphical, intuitive software "dashboard" for easy operation and training.
- Automated GMP reports
 - Complete batch reports
 - Event and logs
 - Chromatograms and calibration history
- Precision Hardware
 - 100:1 turndown ratio for accurate gradients
 - Low system volume
 - Integrated pH, UV, condition flow cell
- Safety features
 - Active air trap minimizes "bubble trap dilution"
 - Leak and pressure alarms
 - Built in safety interlocks
 - Software method checks
 - Software security prevents unauthorized operation and tampering of data

The Right Combination of Simplicity and Power

The simplicity of QuantaSep®'s interactive graphical user interface gives the system its efficiency and power by making it easy to learn and use. You don't need to refer to printed flow diagrams and manuals. Because, even if you have little or no computer experience, QuantaSep®'s intuitive on-screen operations make the system easy to use-from the moment you plug it in..

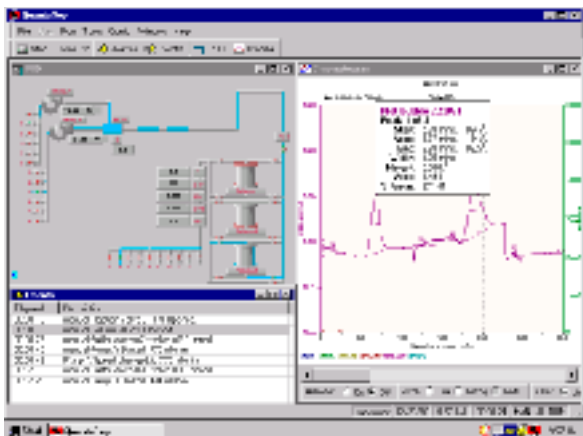

QuantaSep® Has What It Takes

The QuantaSep® System can automate many of the routine manual tasks related to process development and free up key scientists time to handle additional mission-critical projects. It allows the process development departments to manage the deadlines to produce clinical/pre-clinical material, scale-up processes and trouble-shoot or downscale production problems all within existing headcount and budgetary constraints. The payback for these systems has ranged from 3-18 months!



SEPRAGEN
CORPORATION

"INNOVATORS IN SEPARATION TECHNOLOGIES"

1205 San Luis Obispo Avenue, Hayward, CA 94544 U. S. A.
TEL: 510-475-0650 FAX: 510-475-0625
info@sepragen.com HTTP://WWW.SEPRAGEN.COM

QuantaSep® 1000LX Specifications

- 2 x12 user selectable buffers
- 2 x 3 user selectable columns forward/reverse/bypass mode
- 2 x10 user selectable fractions with fraction collect capability
- Sanitary design
- Low flow range: 1-100 ml/min, $\pm 3\%$
High flow range: 10-1000 ml/min, $\pm 3\%$
- Low flow hold up volume: 18 ml
High flow hold up volume: 120 ml
- Small component footprint: 21"X26"X15"
- 5%-95% gradient capability
- Pressure rating: 50 psi
- Low flow tubing: 1/16" I.D.
High flow tubing: 1/4" I.D.

Wetted Materials of Construction

Silicone, polyethylene and polypropylene tubing; polypropylene valves and fittings; PTFE diaphragms; glass/quartz/316L stainless steel sensor components.

Components Specifications

280 nm UV 0-2 Au
254 nm UV 0-2 Au (Optional)
Conductivity Sensor 0-200 ms, $\pm 5\%$
pH Sensor 1-14 pH, ± 0.2
Air Bubble $> 0.20''$
Pressure Sensor 0-100 psi ± 1 psi
Leak Sensor 5 ml leak volume

Warranty, Installation & Support

All QuantaSep®'s are backed by Sepragen's one year limited warranty. Installation and training, along with IQ/OQ are performed by Sepragen Service Group. Preventative maintenance and Total maintenance service is provided at added cost. Validation support packages are also available.

Call (510) 667-1004 to schedule a free confidential ROI/Efficiency assessment of your purification process.